


ZENITH III HyLo ELECTRICAL WIRING

LED FUNCTIONS

LED #	FUNCTION
1	Not Used
2	Post Sw DOWN
3	Not Used
4	KEY Sw
5	Table DOWN Limit
6	Table UP Limit
7	Not Used
8	Not Used
9	Not Used
10	Post Sw UP
11	Ready Sw / E STOP

LOGIC BOARD	POWER BOARD
1	A - DOWN Solenoid Valve
2	E - PUMP Motor
3	
4	
5	
6	
7	
8	
9	
10	
11	

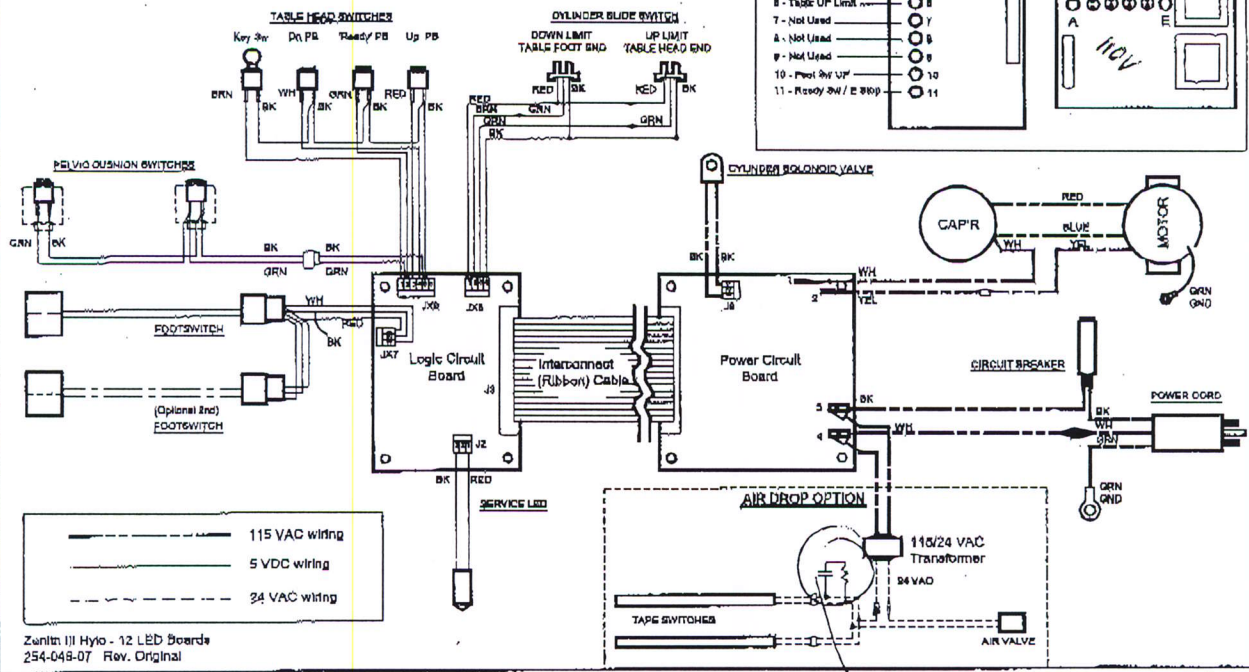

- - - - - 115 VAC wiring
 _____ 5 VDC wiring
 - - - - - 24 VAC wiring

Zenith III HyLo - 12 LED Boards
254-048-07 Rev. Original

E.R.S.