

ZENITH II  
Timer Conversion Kit  
115 V. System

800-537-03  
100-019-00

Panels and Controls


224-031-00

TO

Air-O-Tronics

224-037-01

Material: Timer 224-037-01 (1)  
Connectors 222-004-06 (4)  
Double Faced Tape 252-003-00 (4 inches)  
Insulating Pad 265-002-01 (1)  
(Place on cover over timer)

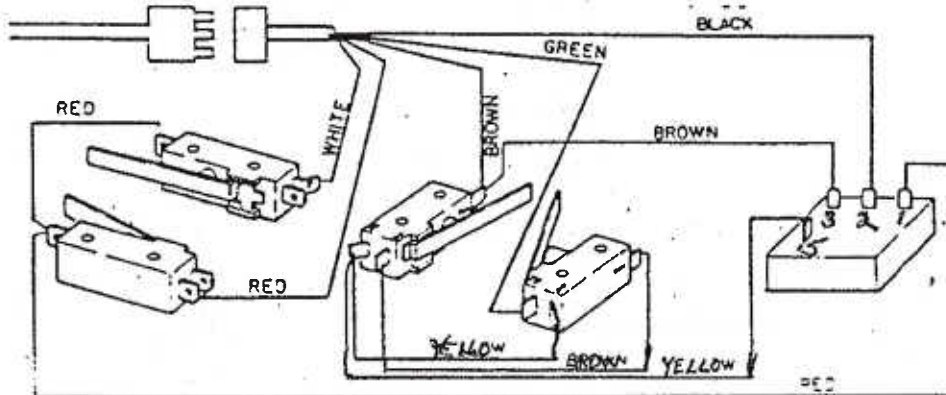

Red # 1

Yellow # 5

Brown # 3

Black # 2

If TABLE still will not RUN  
CHANGE 2 AND 3 AROUND